



Maquet Volista Access II
Surgical Light
Access to more than light

GETINGE 

Maquet Volista Access II Overview

- ⊕ Perform surgery with safety in mind
 - ⊕ Enhance visibility and diagnostic
 - ⊕ Stable illumination
 - ⊕ Prevent burn and tissue desiccation



- ⊕ Improve surgical staff well-being
 - ⊕ Effortless handling
 - ⊕ Peace of mind thanks to dimming sensitivity
 - ⊕ Easy access to lighting controls
 - ⊕ Access to video

Maquet Volista Access II Surgical Light


Perform surgery with safety in mind

Access to adjustable color temperature

Based on a mix of a new generation of white LEDs (warm and cold) Maquet Volista Access II offers the possibility to choose between adjustable color temperature (3,900 ; 4, 200 and 4,500 K) or a fixed color temperature (4,200 K).

The illumination you can rely on



An electronic architecture based on High Efficient Collimators (HECOL) to optimize the useful light throughout the entire surgery. The micro lenses integrated on the collimator homogeneously distribute the lighting inside the patient's cavity.

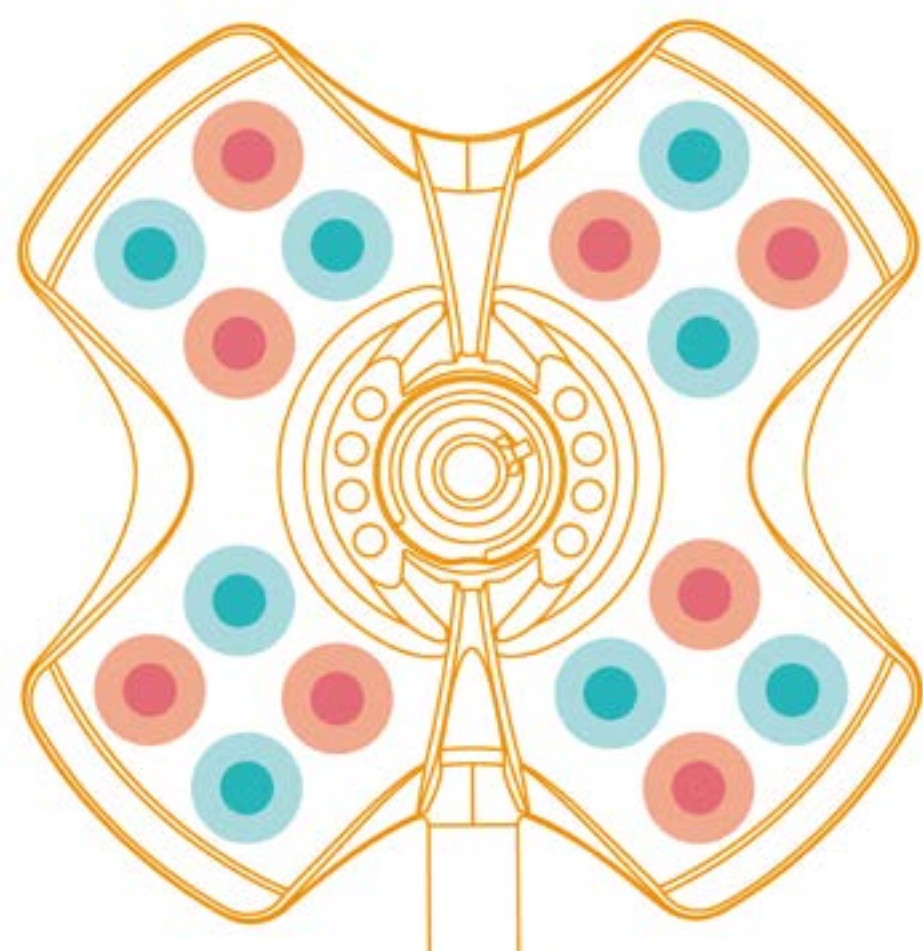

Stable illumination from start to finish

FSP (Flux Stability Program) applies an automatic compensation, cancelling out the heat up and the drop in intensity inherent to LED technology

[Learn more](#)

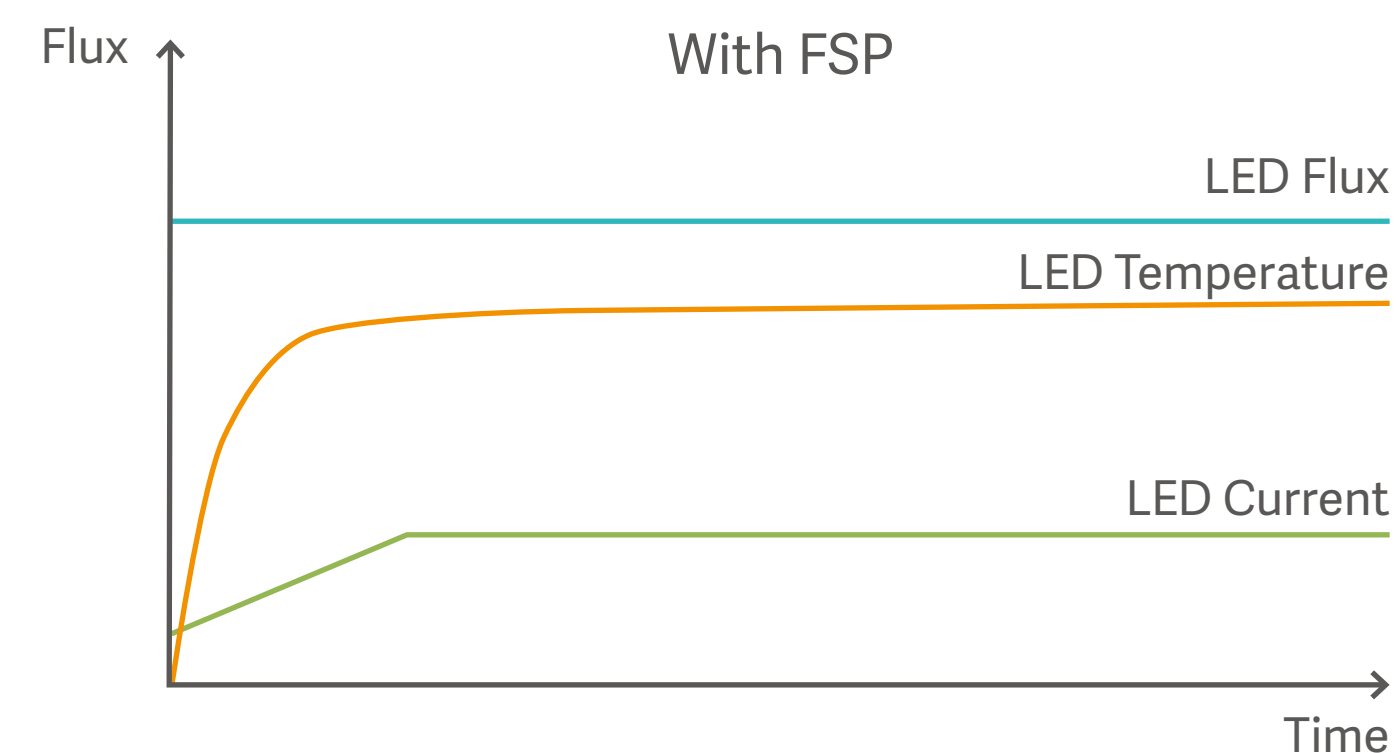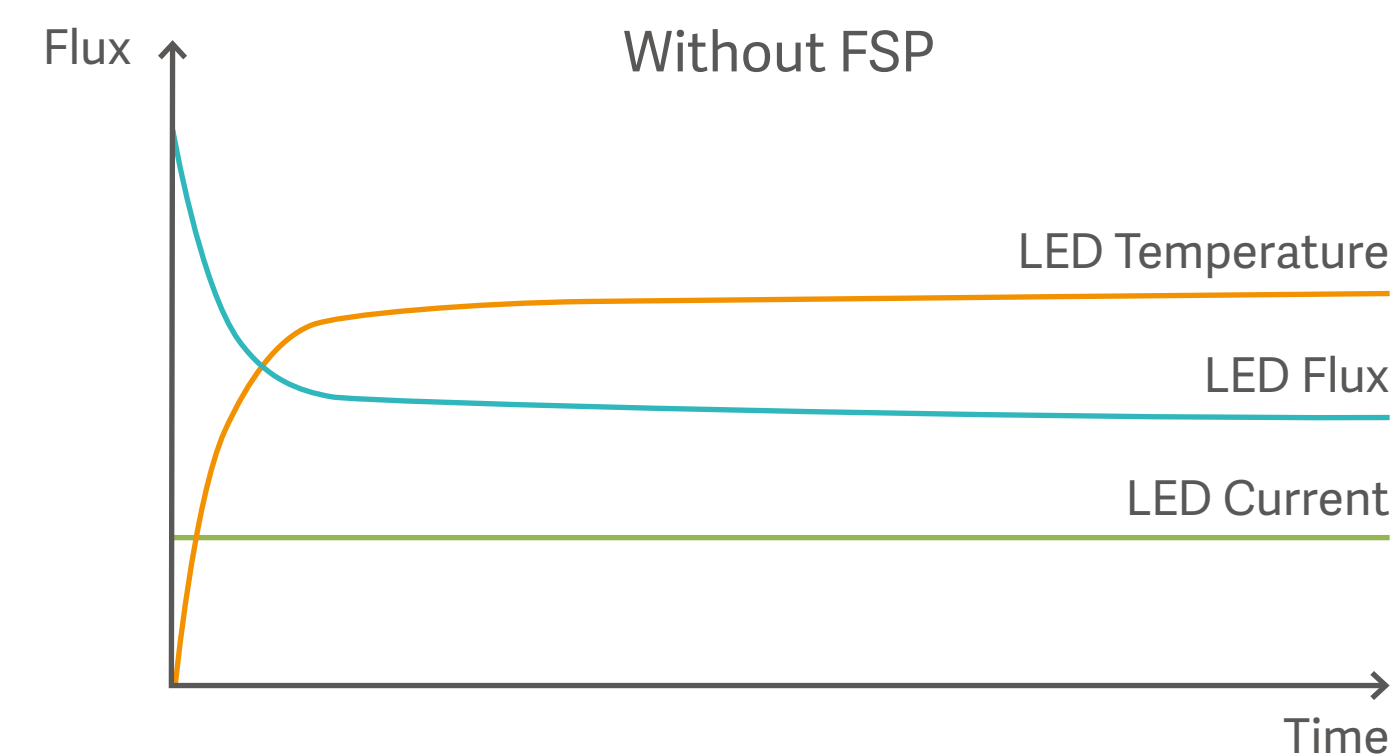

Maquet Volista Access II Surgical Light

Perform surgery with safety in mind

Optimal visibility for minimally invasive surgery


Maquet Volista Access II offers green ambient lighting at the center of the light head to minimize glare and improve contrast on monitors during MIS. The ambient light provides enough illumination to help surgical staff to move safely in the darkened ORs.

Hygienic design

Cleaning made simple thanks to the injected underside and a smooth touch keypad. Maquet Volista has a special paint coating that is designed to reduce bacterial colonization with disinfection.

Safe irradiance

Protect delicate patient tissues thanks to a visual assistance to manage the amount of heat induced by the overlapping light patches in the surgical wound.

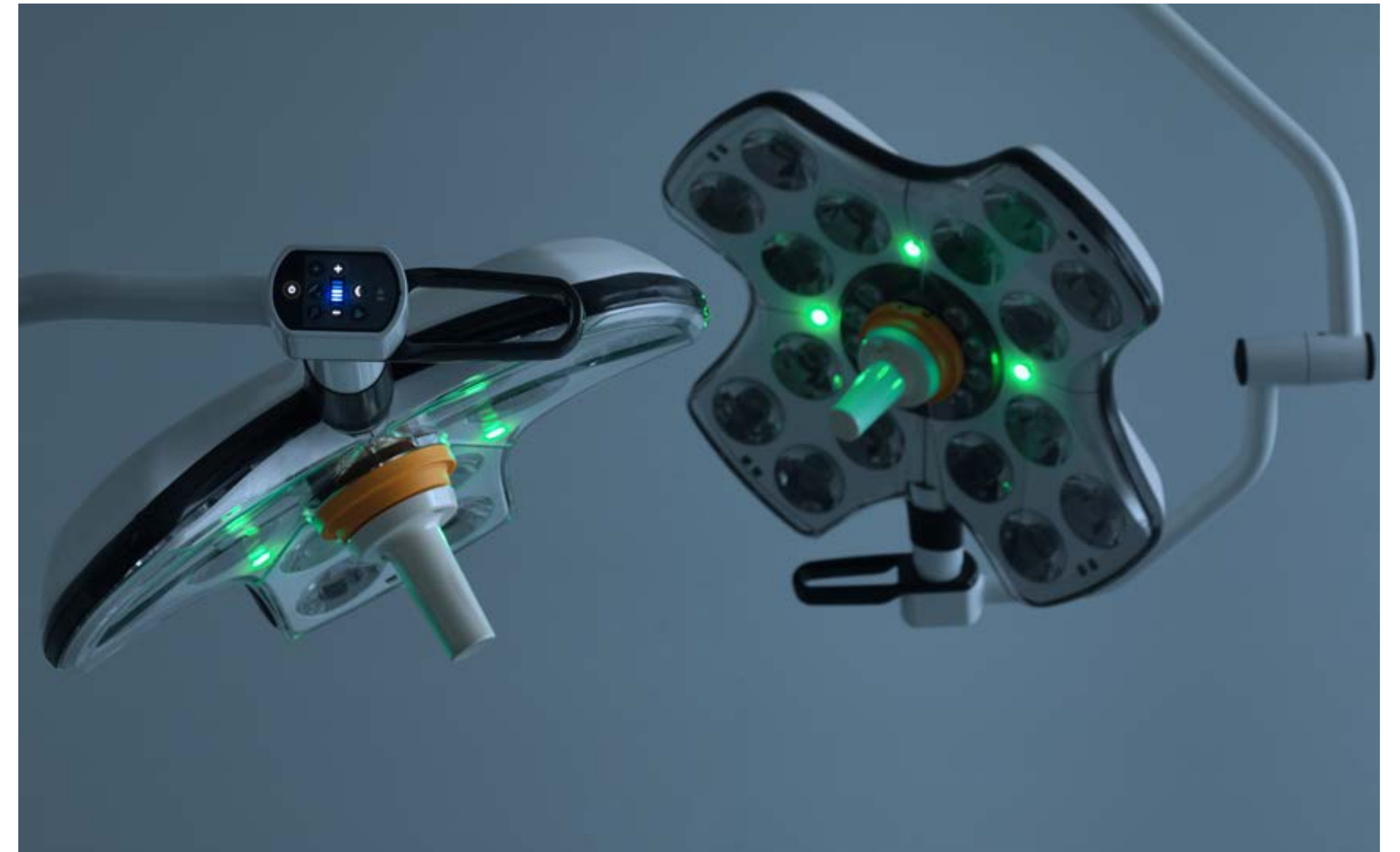

Maquet Volista Access II Surgical Light

Improve surgical staff well-being

Access to video

Full HD wireless cameras can be mounted directly on the light head. The Quick Lock System allows cameras to be quickly and easily connected and disconnected from the light head and to be moved between surgical suites. This tool-free system minimizes the camera's setup time between procedures, and maximizes the use of cameras throughout the surgical suite.

Peace of mind

In the case of a dual configuration, if the color temperature is adjusted on one light head, the second one will be also adjusted due to synchronization.

Dimming sensitivity

With dimming following human eye sensitivity per Fechner's law, each dimming step is seen by human eyes with the same amount of light difference, in order to provide a smooth and adapted range of illumination.



Maquet Volista Access II Surgical Light

Technical data

+ Control panels



Light head keypad




Capacitive wall keypad


+ Optical characteristics

Maximum illumination	160,000 lux
Diameter d10	20-25 cm
Illumination depth at 60%	52 cm
Color temperature	Fixed: 4,200 K Adjustable: 3,900 K / 4,200 K / 4,500 K
Color rendering index	95
FSP (Flux stability program)	Yes

+ Rotation angles

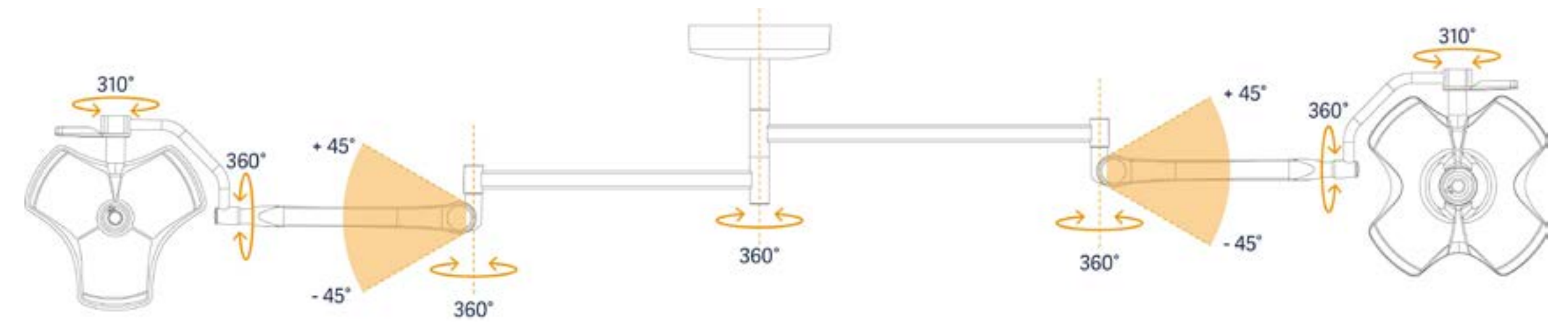

SF version

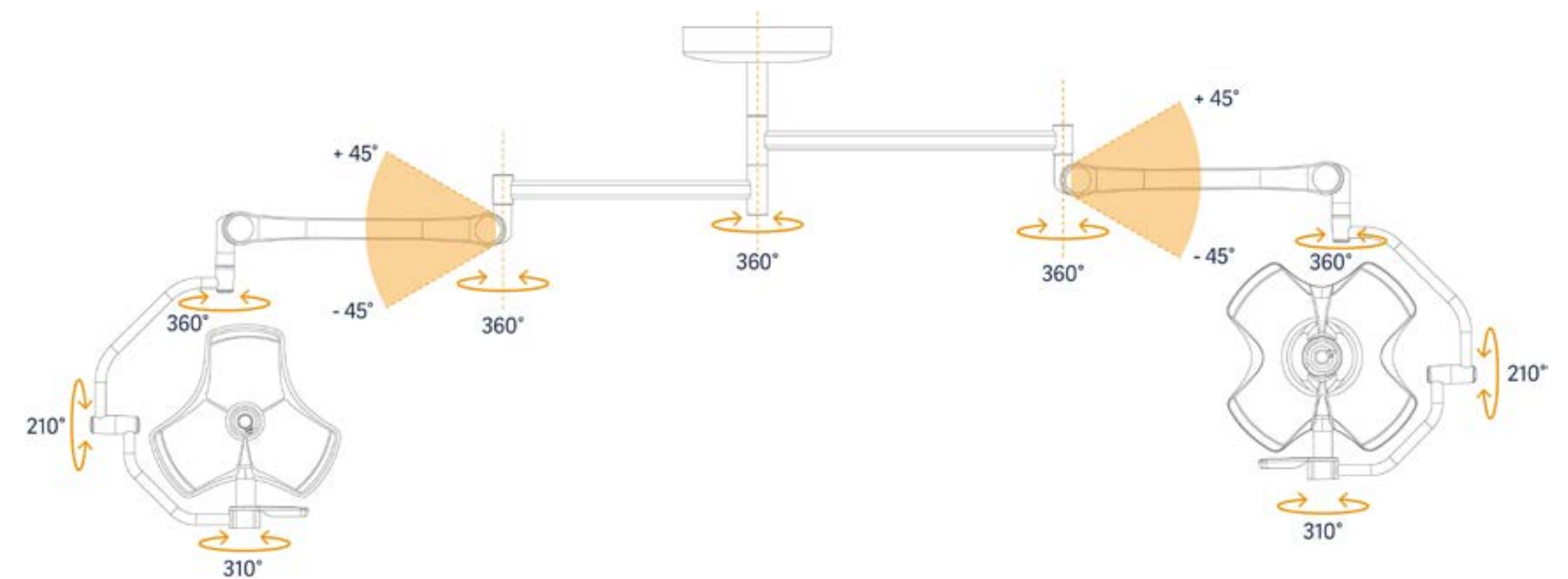


DF version

This information is intended for an international audience outside the US.

This information is aimed exclusively at healthcare professionals or other professional audiences and are for informational purposes only, is not exhaustive and therefore should not be relied upon as a replacement of the Instructions for Use, service manual or medical advice. Getinge shall bear no responsibility or liability for any action or omission of any party based upon this material, and reliance is solely at the user's risk.

Any therapy, solution or product mentioned might not be available or allowed in your country. Information may not be copied or used, in whole or in part, without written permission by Getinge.

Manufacturer: MAQUET SAS · Parc de Limère · Avenue de la Pomme de Pin · CS 10008 Ardon · 45074 Orleans, cedex 2 · France +33 (2) 38 25 88

© 2023 Getinge | Getinge and **GETINGE**  are trademarks or registered trademarks of Getinge AB, its subsidiaries or affiliates. · DMS-0004812 v3 | All rights reserved.

GETINGE 