



Getinge Assured

Comprehensive range of indicators for wash monitoring and sterility assurance



Added assurance

- for clean and sterile goods

Getinge Assured indicators provide you with the added assurance that your equipment is clean and sterile, as expected.

Every day you reprocess a variety of instruments at different times, temperatures and conditions. You need reassurance that your sterile goods have been thoroughly cleaned and sterilized in accordance with your local standards and hospital requirements.

Getinge is a world leader in cleaning, disinfection and sterilization equipment, with solutions for the healthcare and life science sectors. We understand that you need your disinfection and sterilization equipment to consistently perform to manufacturer's specifications in full compliance with all requirements and standards.

Getinge Consumables

- The complete solution

The Getinge Consumables offering is a comprehensive portfolio of products to support you in every step of your workflow. Getinge products conform to relevant standards, ensuring consistent and accurate quality of your results.

Getinge Clean

Getinge Clean is a complete detergent offering for global markets requiring enzymatic, alkaline or neutral cleaning solutions, suitable for both manual and automated processing. We offer specialized chemistry — such as pre-treatment foam to help keep instruments moist. We also offer rinse aid, lubricants, neutralizers and various descaling products including chamber foam as part of our maintenance offering.



Getinge Pack

Keep your patients safe by using high-quality sterilization products. Even the most robust disinfection and sterilization processes can be rendered ineffective if the integrity of a medical pouch or its closure has been compromised during storage, handling, or transportation. That's why Getinge Pack Sterilization Pouches and Rolls have been designed to provide an effective microbial barrier that significantly minimizes the risk of recontamination.



Standard pouch



Standard roll



Pouch made of Tyvek®



Roll made of Tyvek®

DuPont™ and Tyvek® are trademarks or registered trademarks of E.I. du Pont de Nemours and Company.

Get accurate and consistent results

- Every time, every day



Improve patient safety

You can reduce the risk of hospital-acquired infections. Getinge Assured indicators are a relatively low-cost method of monitoring the success or failure of the entire decontamination process. This will help manifest itself in fewer cancelled operations, just-in-time delivery of sterile instruments and contribute to the optimization of the patient's experience.



Increase efficiency in your workflow

Our training programs — available both face-to-face and electronically — give you the tools you need to ensure consistent and accurate results for every point of your workflow. Our sales tools and reference guides make it easy to interpret our products. All products come with clear pass/fail reference guides.



Compliance with relevant international and local standards

Our products give you complete peace of mind. Getinge Assured products are all manufactured in ISO 9001 and ISO 13485 certified Getinge facilities and are tested against all relevant standards such as EN ISO 15883, EN ISO 11140 parts 1 & 4, and EN ISO 11138. We offer certificates of compliance for the complete range, as well as 3rd party independent certificates for most products.



Reduce environmental footprint

Getinge Assured goes beyond wash monitoring and sterility assurance. We also give you complete control in how you impact your local environment. Our unique ink technology is non-toxic and contains no heavy metals in their base ingredients.

Wash Release Monitors

- The first important step in sterile processing

Effective cleaning is the first important step in sterile processing. The best way to verify cleanliness is a reliable testing protocol that includes test soil, wash monitoring, and residual protein testing.

Used routinely in washer-disinfectors or ultrasonic baths, our Getinge Assured Wash Monitors provide independent verification that the cleaning process is consistently providing the desired result.

As a final assurance of cleanliness, we offer a full range of residual protein tests for surfaces and instruments – including hollow lumens and endoscopes. With an instant protein reaction, these tests are quick and easy to use and can detect residual protein down to 1 microgram in less than 10 seconds.



Getinge Assured Wash Monitor



Getinge Assured Wash Monitor F/R



Getinge Assured Wash Monitor Ultrasonic



Getinge Assured Wash Monitor U/L



Getinge Assured Wash Monitor Holder



Getinge Assured Wash Monitor Holder Flexible



Getinge Assured Wash Monitor Holder Rigid



Getinge Assured Test Soil



Getinge Assured Protein Test Instrument Surface



Getinge Assured Protein Test Instrument Lumen



Getinge Assured Protein Test Flexible Endoscope 2,5 m



Getinge Assured Protein Test Dental

Machine Release Monitors

- Assurance of equipment performance

The Bowie-Dick-type test remains the benchmark for daily testing of the mechanical function of your vacuum-assisted steam sterilizer.

Getinge Assured Bowie-Dick tests monitor the effective air removal and steam penetration allowing you to release your sterilizer for daily use. Getinge's innovative Bowie-Dick type tests have been developed and tested to conform to the requirements of EN ISO 11140-4.

With highly sensitive advanced ink, our technology provides greater accuracy and user information on process conditions compared to conventional test packs.



Getinge Assured AdvaTec Bowie-Dick Test Pack 134°C*



Getinge Check Mini Bowie-Dick Pack 134°C**




Getinge Assured Bowie-Dick Test Pack 134°C* & 121°C




Getinge Check Bowie-Dick Test Pack 134°C**



Getinge Assured Bowie-Dick Test Card TT

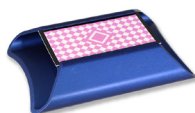

Getinge Assured Bowie-Dick Test Card 134°C

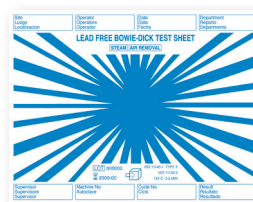

Getinge Assured Bowie-Dick Test Sheet 134°C

* For use in Getinge GSS67H and HS66 Turbo Steam Sterilizers

** For use in all sterilizers except Getinge GSS67H and HS66 Turbo

The compact design of the Getinge Assured AdvaTec and the Getinge Check Mini Bowie-Dick Test requires less paper and minimizes usage cost, transport cost, storage space, and environmental footprint.

Load Release Monitors

- Assurance of effective sterilization conditions

A chemical or biological process challenge device (PCD) gives the assurance that correct conditions have been achieved in the sterilizer chamber.

Product name	Spore Type	Incubation
Superfast 20 Biological Indicator (Steam)	10 ⁶ <i>Geobacillus stearothermophilus</i>	20 minutes
Superfast 20 Biological Indicator (VH2O2)	10 ⁶ <i>Geobacillus stearothermophilus</i>	20 minutes
Smart-Read Biological Indicator (Steam)	10 ⁵ <i>Geobacillus stearothermophilus</i>	10 hrs Fail; Cycle Verification: 3 to 5 hrs
Self-Contained Biological Indicator (Steam)	10 ⁵ <i>Geobacillus stearothermophilus</i>	24 hours
Self-Contained Biological Indicator (EO)	10 ⁶ <i>Bacillus atropheus</i>	48 hours
Self-Contained Biological Indicator (VH2O2)	10 ⁶ <i>Geobacillus stearothermophilus</i>	24 hours



Getinge Assured Superfast 20 Biological Indicator (Steam and VH₂O₂)



Getinge Assured Express Incubator



Getinge Assured Smart-Read Biological Indicator (Steam)



Getinge Assured Smart-Well Incubator with Printer



Getinge Assured Self-Contained Biological Indicator (Steam, EO and VH₂O₂)



Getinge Assured Multi Incubator



Getinge Assured Smart-Read Biological Test Pack



Getinge Assured Self-Contained Biological Test Pack



Getinge Assured Batch Monitoring System 3.5, 4, 5.3 & 7



Getinge Assured Helix Test 3.5, 4, 5.3, 7 & Formaldehyde

Getinge Assured 20 Superfast Biological Indicators for Steam and VH₂O₂ provides users the assurance that the sterilization load is safe to release in just 20 minutes.


Pack Release Monitors

- Sterility assurance at each instrument pack level

When a pack is opened, OR staff can quickly verify sterilization conditions were achieved, and the pack is safe for use.

Getinge Assured Emulating Indicators are the most advanced chemical indicators for monitoring the specific steam sterilization parameters of time, temperature, and steam inside the pack.

Minutes at 134°C	Emulating Indicator Strip	Emulating Indicator Adhesive
3.5		●
4	●	●
5.3	●	●
7	●	●
18	●	

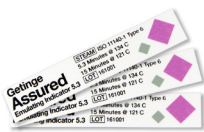

Getinge Assured Emulating Indicators



Getinge Assured Emulating Indicators, Adhesive



Getinge Assured Integrating Indicator



Getinge Assured MI Steam Migrating Integrator Strip




Getinge Assured Multicritical Process Variable Indicator (Steam)



Getinge Assured Multicritical Process Variable Indicator (EO)



Getinge Assured Multicritical Process Variable Indicator (Formaldehyde)

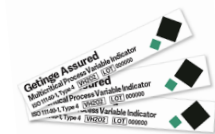

Getinge Assured Multicritical Process Variable Indicator (VH₂O₂)



Getinge Assured Multicritical Process Variable Indicator (Dry Heat)



Getinge Assured Green Tape (Steam) 18mm & 24mm



Getinge Assured Record Labels for Steam, EO & VHP

Label Color	Product name	Use	Order No.
-------------	--------------	-----	-----------

Wash Release

	Test Soil	Washer-disinfectors	504228300
--	-----------	---------------------	-----------

Wash Monitoring and Holders

	Wash Monitor Wash Monitor Wash Monitor	Washer-disinfectors	503864500 503864600 503864700
	Wash Monitor F/R • Flexible & Rigid	Washer-disinfectors	6005500586
	Wash Monitor Holder Wash Monitor Holder R • Rigid Lumen Wash Monitor Holder F • Flexible Lumen	Washer-disinfectors / Ultrasonic	6005500585 503865000 503865100

Wash Monitoring, Ultrasonic

	Wash Monitor U • Ultrasonic Test	Ultrasonic	503865200
	Wash Monitor U/L • Ultrasonic Test	Ultrasonic	6005500587

Protein Tests

	Protein Test Instrument Surface Protein Test Instrument Lumen	Washer-disinfectors / Ultrasonic	503878500 503911200
	Protein Test Flexible Endoscope 2,5 m x 2,8 mm Protein Test Flexible Endoscope 2,5 m x 3,7 mm Protein Test Flexible Endoscope 2,5 m x 1,9 mm Protein Test Dental	Washer-disinfectors / Ultrasonic	503964400 6005500296 6005500297 6005500000

Machine Release

Bowie-Dick Tests

	Getinge Assured Bowie-Dick Test Pack 134°C, Getinge Sterilizers	Steam	504052100
	Getinge Check Bowie-Dick Test Pack 134°C, Sub-trans cycles	Steam	6001154901
	Bowie-Dick Test Pack 121°C	Steam	504052200
	Bowie-Dick Test Card TT (Table Top Sterilizers)	Steam	6005500588
	Bowie-Dick Test Card 134°C	Steam	6005500589
	AdvaTec Bowie-Dick Test Pack, Getinge Sterilizers	Steam	504225100
	Bowie-Dick Test Sheet 134°C	Steam	6005500590
	Getinge Check Bowie-Dick Mini-Pack, Sub-trans	Steam	6001155600

Load Release

Helix Tests

	Helix Test 3.5	Steam	504052700
	Helix Test 4 Helix Test 5.3 Helix Test 7	Steam	504052800 504052900 504053000
	Helix Test Formaldehyde	Formaldehyde	504053100
	Getinge Assured Batch Monitoring System PCD Device	Steam	6005500575

Label Color	Product name	Use	Order No.
	Batch Monitoring System Helix Test 3.5 (Refill Kit)	Steam	6005500037
	Batch Monitoring System Helix Test 4 (Refill Kit)		6005500038
	Batch Monitoring System Helix Test 5.3 (Refill Kit)		6005500039
	Batch Monitoring System Helix Test 7 (Refill Kit)		6005500040
Biological Indicators			
	Superfast 20 Biological Indicator	Steam	6005500506
	Superfast 20 Biological Indicator	VH ₂ O ₂	6005500507
	Express Incubator	20 minutes	6005500511
	Self-Contained Biological Indicator	Steam	504053900
	Self-Contained Biological Indicator	EO	504054000
	Self-Contained Biological Indicator	VH ₂ O ₂	6005500025
	Smart-Read Biological Indicator	Steam	504054100
	Smart-Read Biological Test Pack	Steam	504054500
	Self-Contained Biological Test Pack		504054600
	Smart-Well Incubator with Printer	10 hrs	504054400
	Multi Incubator	24 hrs	504054700

Pack Release

Tape			
	Green Tape 18 mm	Steam	6001154902
	Green Tape 24 mm		6001154903
Record Labels			
	Record Label Steam	Steam	6005500408
	Record Label EO	EO	6005500409
	Record Label VHP	VH ₂ O ₂	6005500410
	Record Label Gun	-	6005500411
Type 4 Indicators			
	Multicritical Process Variable Indicator	EO	504051600
	Multicritical Process Variable Indicator	Steam	504051700
	Multicritical Process Variable Indicator L	Steam	504051800
	Multicritical Process Variable Indicator	Formaldehyde	504051900
	Multicritical Process Variable Indicator	Dry Heat	504054800
	Multicritical Process Variable Indicator	VH ₂ O ₂	6005500509
Type 5 Indicators			
	Integrating Indicator	Steam	6005500594
	MI Steam Migrating Integrator Strip		6005500510
Emulating Indicators - Non adhesive			
	Emulating Indicator 4	Steam	504051000
	Emulating Indicator 5.3		504050800
	Emulating Indicator 7		504050900
	Emulating Indicator 18		504051400
Emulating Indicators - Adhesive			
	Emulating Indicator 3.5 -A	Steam	504051100
	Emulating Indicator 4 - A		504052600
	Emulating Indicator 5.3 -A		504051200
	Emulating Indicator 7 -A		504051300
Miscellaneous.			
	VH ₂ O ₂ Indicator	VH ₂ O ₂	504052000



This document is intended to provide information to an international audience outside of the US.

This information is aimed exclusively at healthcare professionals or other professional audiences and are for informational purposes only, is not exhaustive and therefore should not be relied upon as a replacement of the Instructions for Use, service manual or medical advice. Getinge shall bear no responsibility or liability for any action or omission of any party based upon this material, and reliance is solely at the user's risk.

Any therapy, solution or product mentioned might not be available or allowed in your country. Information may not be copied or used, in whole or in part, without written permission by Getinge.

The products presented may be pending regulatory approvals to be marketed in your country.
Contact your Getinge representative for more information.

Manufacturer:

Getinge Sterility Assurance • 74 Inverness Drive East • Englewood, CO 80112 • United States • +1 303-660-4201 • info@getinge.com

© 2023 Getinge | Getinge and **GETINGE** * are trademarks or registered trademarks of Getinge AB, its subsidiaries or affiliates. | DMS-0002438-rev2 | All rights reserved.

www.getinge.com