

Getinge Pack

Standard Packaging

Self-seal sterilization pouches

Providing high quality and strong microbial barriers for all your processed hospital instruments.

Self-seal sterilization pouches are pre-folded and assure accurate and fast closing without the requirement of a heat sealing device. They are intended to be used to enclose medical devices to be sterilized in steam, ethylene oxide (EO) or formaldehyde sterilizers.


Made of high grade medical paper and multilayer PET/PP copolymer film, Self-seal pouches ensure an effective microbial barrier for all your processed instruments. Water based and non-toxic process indicators printed on the medical grade paper surface provide an accurate visual confirmation by color change and will help to differentiate between processed and unprocessed packages.

A wide range of standard sizes allows you for optimum choice of a correct-sized pack for each instrument and for double packing purposes.

Getinge Pack Self-seal pouches are CE marked and in compliance with international standards ISO 11607-1, ISO 11140-1 and EN 868-5.

This document is intended to provide information to an international audience outside of the US.

www.getinge.com

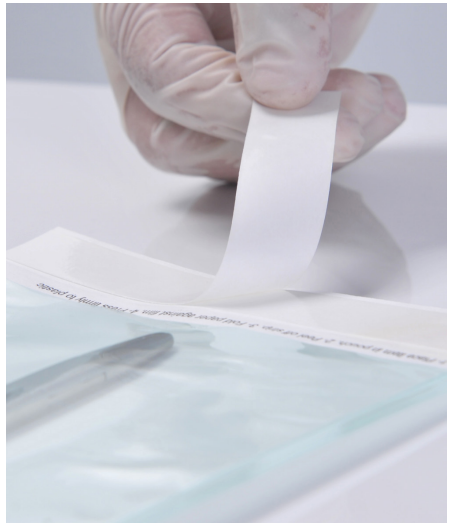

GETINGE 