




OEM Medical Collagen Solutions

from Getinge


OEM Medical Collagen Solutions

Getinge is a global provider of innovative solutions for operating rooms, intensive care units, sterilization departments and for life science companies and institutions since 1904. Based on our firsthand experience and close partnerships with clinical experts & healthcare professionals, Getinge has become a trusted medical OEM supplier worldwide, your partner in the collagen business.

Having a flexible manufacturing process designed to preserve the integrity of the collagen fibers, we are able to partner with you in support of the development of your medical devices.

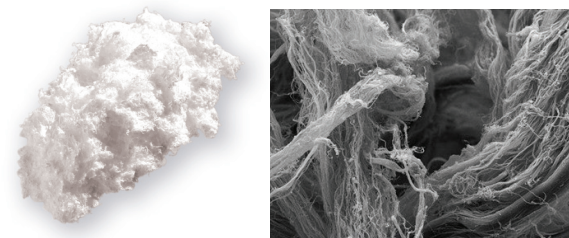

- Micro fibrillar material suitable for collagen slurries
- Effectively used in cavities and irregular surfaces

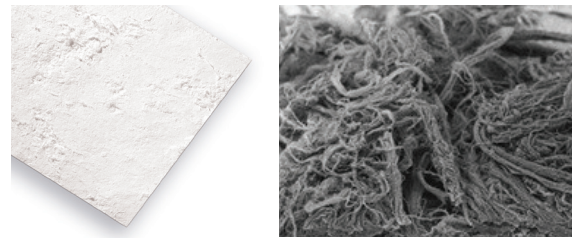

- Compressed/uncompressed sheets with variable density and thickness
- Applicable for small and large wounds

Getinge has been manufacturing Bovine Type I collagen for over 30 years. Our collagen matrix has been used in many clinical applications including vascular closure device and as a coating for textile vascular grafts.



Fibrillar



Uncompressed sheet



Compressed sheet

Collagen source

- Sourced from Achilles tendon of animals younger than 4 years of age
- New Zealand farmed cattle
- GBR-1 rated tendon supplier with certificate of origin from New Zealand Ministry
- Approved supplier and component traceability
- Harvested through a controlled/documentated process

Process environment

- GMP compliant facility
- Class 8 controlled manufacturing clean room facility
- ISO13485-2016
- Quality process designed to preserve the integrity and physical properties of the collagen fibrils
- Flexible manufacturing process to help meet customers' requirements
- Large and small lot capabilities
- Competitive price

Collagen properties

- Wide range of collagen fiber sizes (length, diameter) defines wide spectrum of mechanical properties and resorption profiles
- 100% pure, naturally long collagen fibers combined with short collagen fibrils
- Interwoven strong collagen fibers
- Non toxic type I collagen fibers
- Can remain in the body without risk of reaction
- Adheres to the bleeding surface due to the matrix of long and short fibers stimulating hemostasis

For more information, please contact:

Alain Forget

Director Business Development Collagen

Direct: (201) 995-8861

Mobile: (201) 675-6936

al.forget@getinge.com

Jennifer Herbst

Sr. Global Product Manager Collagen & Cardiac Assist

Direct: (973) 709 7296

Mobile: (201) 574 5684

jennifer.herbst@getinge.com

PN: 0002-08-9800 Rev C · Getinge is a registered trademark of Getinge AB, its subsidiaries, or affiliates in the United States or other countries.

Getinge · 1300 MacArthur Blvd. · Mahwah, NJ 07430 · USA · Phone: +1 201 995 8700 or 1 800 777 4222

www.getinge.com

